


HP-1 α Rabbit Polyclonal Antibody(HP-1 a)

Catalog No: RA20278


Basic Information	
Host species	Rabbit
Applications	WB,IHC
Species Cross-Reactivity	H,R,M
Specificity	Antibody can detects endogenous HP-1αproteins.
Recommended dilutions	WB: 1:1,000-2,000 IHC: 1:200-500
	Optimal dilutions should be determined by the end user.
Applications	
Formulation	Antigen Affinity Purified IgG in PBS, pH 7.4, containing 0.02% sodium azide as Preservative and 50% Glycerol.
Storage	Store at -20°C. Avoid repeated freeze-thaw cycles.
Concentration	1 mg/ml
Clonality	Polyclonal
Background	
Alternative Names	Antigen p25, CBX5, HEL25, Heterochromatin protein 1
Observed band	22
Human Gene ID	23468
Human Swiss-Prot Number	P45973
Background	Heterochromatin protein 1 (HP1) is a family of heterochromatic adaptor molecules involved in both gene silencing and higher order chromatin structure. All three HP1 family members (α , β , and γ) are primarily associated with centromeric heterochromatin.



Selected Validation Data

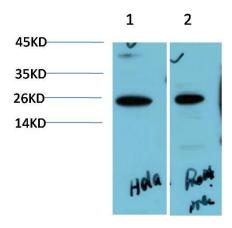

Western blot analysis of 1) Hela Cell Lysate, 2) Rat Brain Tissue Lysate using HP-1a Rabbit pAb diluted at 1:2000.

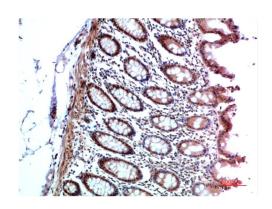

Immunohistochemical analysis of paraffin-embedded Human Colon Carcinoma Tissue using HP-1 α Rabbit pAb diluted at 1:500.