


PCNA Mouse Monoclonal Antibody(12D10)

Catalog No: RCA23


Basic Information	
Host species	Mouse
Applications	WB, IHC
Species Cross-Reactivity	H, R, M
Specificity	Antibody can detects endogenous PCNA protein.
Recommended dilutions	WB: 1:5,000 IHC: 1:200
	Optimal dilutions should be determined by the end user.
Applications	
Formulation	PBS, pH 7.4, containing 0.02% sodium azide as Preservative and 50% Glycerol.
Storage	Store at -20°C. Avoid repeated freeze-thaw cycles.
Concentration	1 mg/ml
Clonality	Monoclonal
Background	
Alternative Names	
Observed band	29
Human Gene ID	5111
Human Swiss-Prot Number	P12004
	PCNA is a marker for cells in early G1 phase and S phase of the cell cycle. It is
	found in the nucleus and is a cofactor of DNA polymerase delta. PCNA acts as a
	homotrimer and helps increase the processivity of leading strand synthesis during
Background	DNA replication. In response to DNA damage, PCNA is ubiquitinated and is
	involved in the RAD6 dependent DNA repair pathway. Two transcript variants
	encoding the same protein have been found for PCNA. Pseudogenes of this gene
	have been described on chromosome 4 and on the X chromosome.



Selected Validation Data

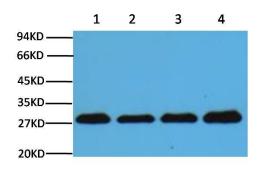

Western blot analysis of Hela (1), Rat brain (2), NIH 3T3 (3), 293T(4) with PCNA mouse mAb(12D10) diluted at 1:5000.

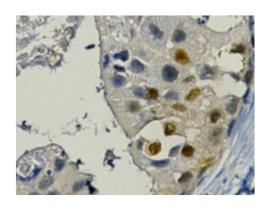

IHC staining of Human breast cancer tissue paraffin-embedded with PCNA mouse mAb (12D10) diluted at 1:200