

α -tubulin Mouse Monoclonal Antibody (Zebrafish Specific)


Catalog No: RCA34

Basic Information	
Host species	Mouse
Applications	WB
Species Cross-Reactivity	Zebrafish
Specificity	The α -tubulin antibody can detects Zebrafish endogenous α -tubulin protein.
Recommended dilutions	WB: 1:5,000-10,000
	Optimal dilutions should be determined by the end user.
Applications	
Formulation	PBS, pH 7.4, containing 0.02% sodium azide as Preservative and 50% Glycerol.
Storage	Store at -20°C. Avoid repeated freeze-thaw cycles.
Concentration	1 mg/ml
Clonality	Monoclonal
Background	
Alternative Names	alpha tubulin 2, B alpha 1, H2 alpha, Hum a tub 1, TUBA1, TUBA1A,Tubulin alpha
	1chain
Observed band	52
Human Gene ID	7277
Human Swiss-Prot Number	P68363
Background	Tubulin is one of several members of a small family of globular proteins. The tubulin superfamily includes five distinct families, the alpha-, beta-, gamma-, delta-, and epsilon-tubulins. The most common members of the tubulin family are α -tubulin and β -tubulin, the proteins that make up microtubules. Each has a molecular weight of approximately 55 KD. Microtubules are assembled from dimers of α - and β -tubulin.



Antibodies · ELISA kits · Biochemical reagents · Recombinant proteins · Immune related reagents

Selected Validation Data

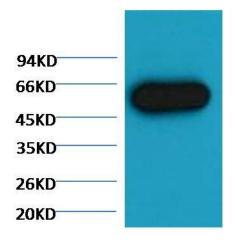

Western blot analysis of Zebrafish skeletal muscle with α tubulin Mouse mAb(Zebrafish Specific) diluted at 1:5,000.