


Antibody(Zebrafish Specific)

Catalog No: RCA35

Basic Information	
Host species	Mouse
Applications	WB
Species Cross-Reactivity	Zebrafish
Specificity	The β I tubulin antibody can recognize Zebrafish endgenousβ I tubulin proteins.
Recommended dilutions	WB: 1:5,000-10,000
	Optimal dilutions should be determined by the end user.
Applications	
Formulation	PBS, pH 7.4, containing 0.02% sodium azide as Preservative and 50% Glycerol.
Storage	Store at -20°C. Avoid repeated freeze-thaw cycles.
Concentration	1 mg/ml
Clonality	Monoclonal
Background	
Alternative Names	beta I tubulin, TBB1, TUBB1, Tubulin beta 1, Class VI beta tubulin
Observed band	50
Human Gene ID	81027
Human Swiss-Prot Number	Q9H4B7
Background	Tubulin is one of several members of a small family of globular proteins. The tubulin superfamily includes five distinct families, the alpha-, beta-, gamma-, delta-, and epsilon-tubulins. The most common members of the tubulin family are α -tubulin and β -tubulin, the proteins that make up microtubules. Each has a molecular weight of approximately 55 KD. Microtubules are assembled from dimers of α - and β -tubulin.



Selected Validation Data

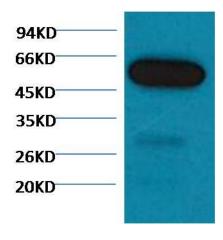

Western blot analysis of Zebrafish skeletal muscle with β I tubulin Mouse Monoclonal Antibody(Zebrafish Specific) diluted at 1:5,000.