

β -actin Mouse Monoclonal Antibody (Zebrafish Specific)


Catalog No: RCA36

Basic Information	
Host species	Mouse
Applications	WB, IHC, IF
Species Cross-Reactivity	Zebrafish
Specificity	Antibody can detects Zebrafish endogenous β-actin protein.
Recommended dilutions	WB: 1:2,000-5,000 IHC: 1:200 IF: 1:200 Optimal dilutions should be determined by the end user.
Applications	
Formulation	PBS, pH 7.4, containing 0.02% sodium azide as Preservative and 50% Glycerol.
Storage	Store at -20°C. Avoid repeated freeze-thaw cycles.
Concentration	1 mg/ml
Clonality	Monoclonal
Background	
Alternative Names	actin antibody, b-actin antibody, beta-actin antibody
Observed band	43
Human Gene ID	60
Human Swiss-Prot Number	P60709
Background	β -actin is one of six different actin isoforms that have been identified. The actin molecules found in cells of various species and tissues tend to be very similar in their immunological and physical properties. Therefore, antibodies against β -actin are useful as loading controls for Western Blotting. However it should be noted that levels of β -actin may not be stable in certain cells. For example, expression of β -actin in adipose tissue is very low and therefore β -actin should not be used as loading control for these tissues.



Antibodies · ELISA kits · Biochemical reagents · Recombinant proteins · Immune related reagents

Selected Validation Data

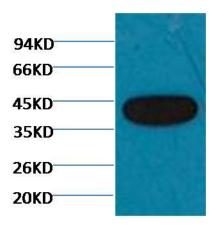

Western blot analysis of Zebrafish skeletal muscle with β -actin Mouse Monoclonal Antibody (Zebrafish Specific) diluted at 1:3,000.

www.reedbiotech.com