


Drosophila Melanogaster GAPDH Rabbit Polyclonal Antibody

Catalog No: RCA69

Dania Information	
Basic Information	
Host species	Rabbit
Applications	WB
Species Cross-Reactivity	Drosophila Melanogaster (H, M, R)
Specificity	Antibody can detects endogenous Drosophila Melanogaster GAPDH protein, predict Cross-Reactivity with Human, Mouse, Rat.
Recommended dilutions	WB:1:5,000-10,000
	Optimal dilutions should be determined by the end user.
Applications	
Formulation	Antigen Affinity Purified IgG in PBS, pH 7.4, containing 0.02% sodium azide as Preservativeand 50% Glycerol.
Storage	Store at -20°C. Avoid repeated freeze-thaw cycles.
Concentration	1 mg/ml
Clonality	Polyclonal
Background	
Alternative Names	G3PDH antibody, BARS-38 antibody, Glycerasldehyde 3 phosphate dehydrogenase antibody
Observed band	38
Human Gene ID	35728
Human Swiss-Prot Number	P07486
	Glyceraldehyde 3 phosphate dehydrogenase (GAPDH) is well known as one of the key enzymes involved in glycolysis. GAPDH is constitutively expressed in almost all
Background	tissues at high levels, therefore antibodies against GAPDH are useful as loading controls for Western Blotting. Some physiological factors, such as hypoxia and diabetes, increase GAPDH expression in certain cell types.



Selected Validation Data

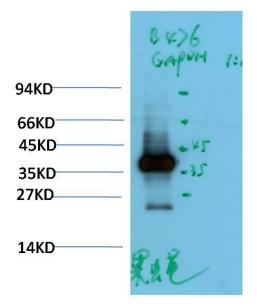

Western blot analysis of Drosophila Melanogaster whole body Lysate using Drosophila Melanogaster GAPDH Rabbit pAb diluted at 1:10,000.